



CATCH PLUS

REVASCULARISATION RANGE



L'INNOVATION SUR MESURE*
www.balt.fr

* innovation, made to measure.

CATCH PLUS

REVASCULARISATION RANGE

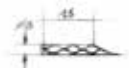
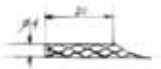
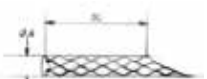

NEW CATCH+ PORTFOLIO

- Complete range of revascularisation thrombectomy devices for the treatment of Acute Ischemic Stroke.
- The "Stentriever" design allows catching the clot within the cells and restore the blood flow.

NEW DESIGN



- Self-expandable laser-cut stent made of nitinol.
- Consistent radial force.
- Mounted on its own delivery wire (length: 180 cm).
- 4 radio-opaque markers : 3 distal (A) + 1 proximal (B).

| COMMERCIAL NAME | REFERENCE | STENT | | CHARACTERISTICS | UNCONSTRAINED STENT | | COMPATIBLE DELIVERY MICROCATETER |
|---|-----------------------|----------------|-------------|---|---------------------|-------------|--|
| | | NOMINAL Ø (MM) | LENGTH (MM) | | Ø (MM) | LENGTH (MM) | |
| CATCH+ MINI  | CATCHMINI 3X15 | 3 | 15 | MINI: for vessels down to Ø1,5 mm | 4 | 15 | VASCO+10 or any 10/14 microcatheter with ID min .017" |
| CATCH+  | CATCH | 4 | 20 | STANDARD: for vessels up to 3 mm | 5 | 20 | VASCO+18 or any microcatheter with ID min .021" |
| CATCH+ MAXI  | CATCHMAXI 6X30 | 6 | 30 | MAXI: for vessels up to 6 mm | 6.5 | 30 | VASCO+21 or any microcatheter with ID min .024" |
| CATCH+ MEGA  | CATCH9 | 9 | 55 | MEGA: for vessels up to 9 mm | 11 | 55 | FARGOMINI or VASCO +35, or any microcatheter with ID min .040" |



10, rue de la Croix Vigneron
95160 Montmorency • France
Tél. : +33 (0)1 39 89 46 41
Fax : +33 (0)1 34 17 03 46
www.balt.fr

NB: All the CATCH+ references are compatible with the use of the balloon tipped guiding catheter CORAIL8F+ ASPI, and also with VASCO+35 ASPI, FARGO+6F ASPI, FARGOMINI.

Delivery microcatheter, aspiration microcatheter and aspiration guiding catheter are not provided with the CATCH+ devices.